


Barron Radial Vacuum Trephine for the recipient cornea



CAT. NO.	DIAMETER
K20-2050	6.0mm
K20-2052	6.5mm
K20-2054	7.0mm
K20-2055	7.25mm
K20-2056	7.5mm
K20-2057	7.75mm
K20-2058	8.0mm
K20-2059	8.25mm
K20-2060	8.5mm
K20-2061	8.75mm
K20-2062	9.0mm

sterile, disposable

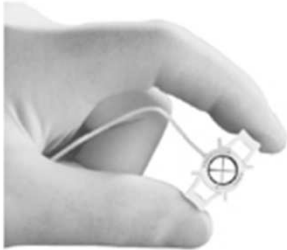


Each sterile, disposable trephine is supplied completely assembled and ready for use.



The Barron Radial Vacuum Trephine is designed for cutting the recipient cornea. It features a 360° suction chamber, divided into 16 compartments. By dividing the suction chamber,

less corneal tissue is aspirated, thus reducing corneal distortion. The result is a consistently accurate, uniform and perpendicular cut. The 16 radial divisions leave a symmetrical pattern on the recipient bed for accurate alignment of sutures to minimize surgically induced astigmatism. These divisions may be inked with a sterile marking pen to be clearly visible for the duration of the surgical procedure. Fine cross hairs are provided for accurate alignment with the visual axis of the recipient. A precisely calibrated rotating mechanism allows the surgeon to advance the blade depth 0.25mm with each 360° rotation. Thus, the trephine is ideally suited for lamellar or penetrating keratoplasty.

RIVENDITORE ESCLUSIVO PER LOMBARDIA E PIEMONTE

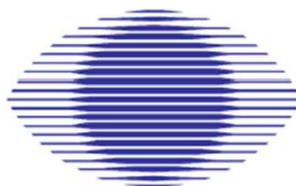

SIR oftalmica

Strumenti chirurgici per oftalmologia
Viale Rosselli nr. 31 — 22100 COMO
Tel. 031-570869 fax 031-572826
E-mail: info@siroftalmica.com

Barron Non-Vacuum Punch for the donor cornea

sterile, disposable

CAT. NO.	DIAMETER	CAT. NO.	DIAMETER
K20-2070	6.0mm	K20-2078	8.0mm
K20-2072	6.5mm	K20-2079	8.25mm
K20-2073	6.75mm	K20-2080	8.5mm
K20-2074	7.0mm	K20-2081	8.75mm
K20-2075	7.25mm	K20-2082	9.0mm
K20-2076	7.5mm	K20-2084	9.5mm
K20-2077	7.75mm		



CUTTING
BLOCK TREPHINE
BLADE




Barron

Artificial Anterior Chamber
holds a (14-18mm) donor corneal graft
sterile disposable


To cut the donor cornea from the epithelial side.

RIVENDITORE ESCLUSIVO PER LOMBARDIA E PIEMONTE

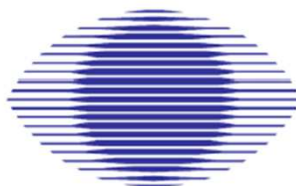

SIR oftalmica

Strumenti chirurgici per oftalmologia
Viale Rosselli nr. 31 — 22100 COMO
Tel. 031-570869 fax 031-572826
E-mail: info@siroftalmica.com

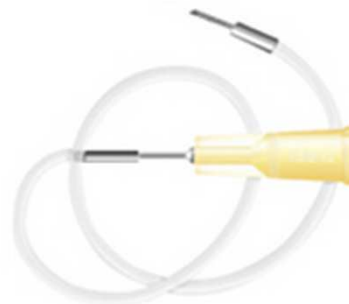


Chamber Maintainer

threaded tip
sterile, disposable (bx/10)

K20-3271	25 gauge
K20-3273	23 gauge
K20-3276	20 gauge

RIVENDITORE ESCLUSIVO PER LOMBARDIA E PIEMONTE

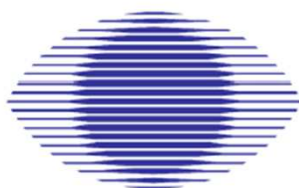

SIR oftalmica

Strumenti chirurgici per oftalmologia
Viale Rosselli nr. 31 — 22100 COMO
Tel. 031-570869 fax 031-572826
E-mail: info@siroftalmica.com